



Recombinant Fc Fragment Of IgG Binding Protein (FcγBP)

Catalog No.: TP10024 50μg

Sequence Information

Species: Human

Gene ID:8857

Swiss Prot:Q9Y6R7

Synonyms:FC(GAMMA)BP

Residues:Cys5176-Ala5344

CTHNGRYLPVNSSLLTSDCSERCSRSSSGLTCQAAGCPPGRVCEVKAECARNCWATRGLCVLSVGANLTTFDGARGATTSPGVYELSSRCPGLQ
NTIPWYRVVAEVQICHGKTEAVGQVHIFFQDGMVTLTPNKGWVWNGLRVDLPAEKLASVSVSRTP

Product Information

Source: Recombinant expression.

Host: *E. coli*

Tags: N-terminal His-Tag

Subcellular Location: Extracellular matrix.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 5% Trehalose.

Original Concentration: 200μg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.3

Predicted Molecular Mass: 19.6kDa

Accurate Molecular Mass: 20kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

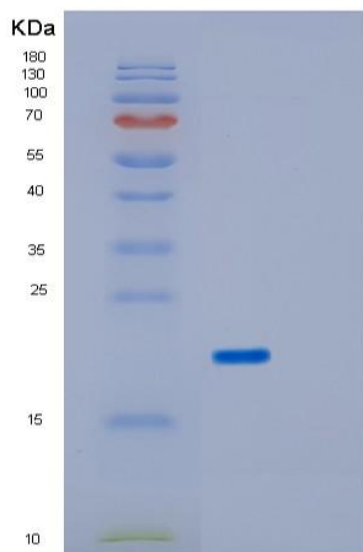


Figure . SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.