

Recombinant Human Matrix Metalloproteinase 1 (MMP1)

Catalog No.: TP09971

50µg

Sequence Information

Species: Human

Gene ID:4312

Swiss Prot:P03956

Synonyms:interstitial collagenase isoform 1
preproprotein , CLG, CLGN

Residues:Phe20-Asn469

FPATLETQEQDVLVQKYLEKYNNLKNDGRQVEKRRNSGPVVEK LKQM QE FGL
KVTGKPDAETLKVMKQPRCGVPDVAQFVLTEGNPRWEQTHLTYRIENYTPDLPR
ADV DHAIEKAFQLWSNV TPLTFTKVSEGGADIMISFVRGDHRD NSPFDGPGGNL
AHAFQPGPGIGGDAHFDEDERWTNNFREYNLHRVAAHELGHSLGLSHSTDIGAL
MYPSTYFSGDVQLAQDDIDGIQAIYGRSQNPVQPIGPQTPKACDSKLTFDAITT
IRGEVMFFKDRFYMRNPFYPEVELNFI SVFWPQLPNGLEAAYEFADRDEV RFF
KGNKYWAVQGQNVLHGYPKDIYSSFGFPRTVKHIDAALSEENTGKTYFFVANKY
WRYDEYKRSM DPGYPKMIAHDFPGIGHKVD AVFMKDGFFYFFHGT RQYKFDPKT
KRILTLQKANSWFNCRKN

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Secreted, Extracellular matrix

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.9

Predicted Molecular Mass: 58.2kDa

Accurate Molecular Mass: 58kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

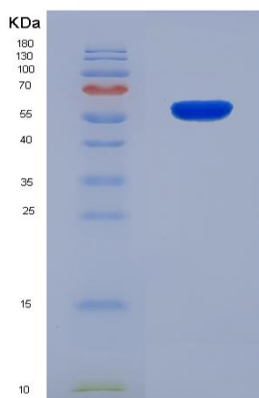


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.