

## Recombinant Osteonectin (ON)

Catalog No.: TP10029

100µg

### Sequence Information

**Species:** Human

**Gene ID:**6678

**Swiss Prot:**P09486

**Synonyms:**SPARC; BM-40; Secreted

Protein,Acidic,Cysteine-Rich;

Basement-membrane protein 40;

Secreted protein acidic and rich in  
cysteine

**Residues:**Ala18~Ile303

APQQEALPDETEVVEETVAEVTEVSVGANPVQVEVGEFDDGAEETEEVVAENP  
CQNHCKKHGKVCELDENTPMVCVCQDPTSCPAIGFEKVCSDNKNKTFDSSCHF  
FATKCTLEGTKKGHKLHLDYIGPKYIPCLDSELTEFPLMRDWLKNVLVTLY  
ERDEDNNLLTEKQKLRVKKIHENEKRL EAGDHPVELLARDFEKNYNMYIFPVHW  
QFGQLDQHPIDGYLSHTELAPLRAPLIPMEHCTTRFFETCDLNDKYIALDEWA  
GCFGIKQKDIDKDLVI

### Product Information

**Source:** Recombinant expression.

**Host:** *E.coli*

**Tags:** N-terminal His-Tag

**Subcellular Location:** Membrane, Secreted, Extracellular matrix.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.0

**Predicted Molecular Mass:** 36.4kDa

**Accurate Molecular Mass:** 37kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

## **[ STORAGE AND STABILITY ]**


**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## **[ IDENTIFICATION ]**

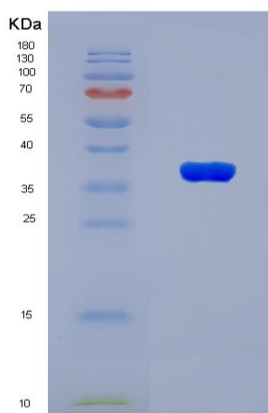

**Figure 1. SDS-PAGE**

## **[ IMPORTANT NOTE ]**

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.