

Recombinant Ubiquitin-like modifier-activating enzyme 1 (UBA1)

Catalog No.: TP11020

100µg

Sequence Information

Species: Human

Gene ID:7317

Swiss Prot:P22314

Synonyms:

Residues:Met1-Arg1058

MSSSPLSKRRVSGPDPKPGSNCSPAQSVLSEVPSVPTNGMAKNGSEADIDEG
YSRQLYVLGHEAMKRLQTSSVLVSGRLGLGVEIAKNILGGVKAVTLHDQGT
WADLSSQFYLR EEDIGKNRAEVSQPR LAELNSYVPVTAYTGPLVEDFLSGFQV
VLTNTPLEDQLRVGEFCHNRGIKLVVADTRGLFGQLFCDFGEEMILTDSNGEQ
LSAMVSMVTKDNPGVVTCLDEARHGFESGDFVSFSEVQGMVELNGNQPM
EIKVLGPYTF SICDTSNFSDIYRGGIVSQVKVPKKISFKSLVASLAEPDFV
VTDFAKFSRPAQLHIGFQALHQFCAQHGRPPRPRNEEDAAELVALAQAVNARALPAVQ
QNNLDEDLIRKLAYVAAGDLAPINAFIGGLAAQEVMKACSGKFMPI
MQWLYFDALECLPEDKEVLTEDKCLQRQNRDYGQVAVFGSDLQEKLGKQKYFLV
GAGAIGCELLKNFAMIGLGC GEGGEIIVTDMDTIEKSNLNRQFLFRP
WDVTKLKSDTAAAVRQMNPHIRVTS HQNRVGP DTERIYDDDFQNL
DGVANALDNVDARMYMDRRCVYYRKP LLESGTLGTKGNVQV
VIPFLTESYSSSQDPPEKSIPICTLKNFPNAIEHTLQWARDEF
EGLFKQPAENVNQYLTDPKFVERTLRLAGTQPLEVLEAVQRS
LVLQRPQTWADCVTWACHHWTQYSNNIRQLLHNFPDPQLTSSGAP
FWSGPKRCPHPLTFDVNNPLHLDYVMAAANLFAQTYGLTGSQDRAA
VATFLQSVQVPEFTPKSGVKIHVS DQELQSANASVDDSRLEELKAT
LPSPDKLPGFKMYPIDFEKDDDSNFHMDFIVAASNLRAENYDIPSAD
RHKSKLIAGKIIPAIAITTTAAVVGVLVCLLEYKVVQGHRQLDSYKNG
FLNALPFFGFSEPLAAPRHQYYNQEWTLWDRFEVQGLQPNGEEMTL
KQFLDYFKTEHKLEITMLSQGVSMLYSFFMPAAKKERLDQPMTEIVSRV
SKRKLGRHVRALVLELCCNDESGEDVEVPYVRYTIR

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.5

Predicted Molecular Mass: 124.3kDa

Accurate Molecular Mass: 124kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

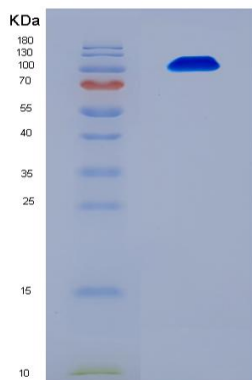


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any



issue if the kit was used in clinical diagnostic or any other procedures.