

Recombinant Suprabasin (SBSN)

Catalog No.: TP11041

100µg

Sequence Information

Species: Human

Gene ID:374897

Swiss Prot:Q6UWP8

Synonyms:

Residues:Asp26-Pro590

DDPIEKVIEGINRGLSNAEREVGKALDGINSGITHAGREVEKVFNGLSNMGSHT
GKELDKGVQGLNHGMDKVAHEINHIGIQAGKEAEKLGHGVNNAAGQVGKEADKL
IHHGVHHGANQAGSEAGKFGQGVDNAAGQAGNEAGRFGQGVHHAAGQAGNEAGR
FGQGVHHAAGQAGNEAGRFGQGAHHGLSEGWKETEKFGQGIHHAAGQVGKEAEK
FGQGAHHAAGQAGNEAGRFGQGVHHGLSEGWKETEKFGQGVHHTAGQVGKEAEK
FGQGAHHAAGQAGNEAGRFGQGAHHAAGQAGNEAGRFGQGVHHGLSEGWKETEK
FGQGVHHAASQFGKETEKLGHGVHHGVNEAWKEAEKFGQGVHHAASQVGKEEDR
VVQGLHHGVSQAGREAGQFGHDIHHTAGQAGKEGDI AVHGVQPGVHEAGKEAGQ
FGQGVHHTLEQAGKEADKAVQGFHTGVHQAGKEAEKLGQGVNHAADQAGKEVEK
LGQGAHHAAGQAGKELQNAHNGVNQASKEANQLLNGNHQSGSSSHQGGATTTPPL
ASGASVNTPFINLPALWRSVANIMP

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.0

Predicted Molecular Mass: 64.4kDa

Accurate Molecular Mass: 64kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

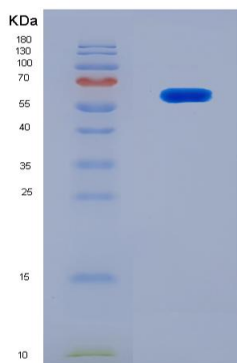


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.