

Recombinant Exportin-1(XPO1)

Catalog No.: TP11077

100µg

Sequence Information

Species: YEAST

Gene ID:853133

Swiss Prot:P30822

Synonyms: CRM1

Residues:Met1-Asp1084

MEGILDFSNDLDIALLDQVVSTFYQGSGVQQKQAQEILTKFQDNPDWQKADQI
LQFSTNPQSKFIALSILDKLITRKWLLPNDHRIGIRNFVVGMIISMCQDDEVF
KTQKNLINKSDLTLVQILKQEWPNWPEFIPELIGSSSSSVNVCENNMIVLKLL
SEEVDFDSAEQMTQAKALHLKNSMSKEFEQIFKLCFQVLEQGSSSSLIVATLES
LLRYLHWIPYRYIYETNILELLSTKFMTSPDTRAITLKCLTEVSNLKIPQDNDL
IKRQTVLFFQNTLQQIATSVMPVTADLKATYANANGNDQSFLQDLAMFLTTYLA
RNRALLESDLESRELLLNHAHQYLIQLSKIEERELFKTTLDYWHNLVADLFYEVQ
RLPATEMSPLIQLSVGSQAISTGSGALNPEYMKRFPLKKHIYEEICSQLRLVII
ENMVRPEEVLVVENDEGEIVREFVKESDTIQLYKSEREVLVYLTHLNVIDTEEI
MISKLARQIDGSEWSWHNINTLSWAIGSISGTMSSEDTEKRFVVTVIKDLLDLTV
KKRGKDNKAVVASDIMYVVGQYPRFLKAHWNFLRTVILKLFEFMHETHEGVQDM
ACDTFIKIVQCKYHFVIQQPRESEPFIQTIIRDIQKTTADLQPQQVHTFYKAC
GIIISEERSVAERNRLSDLMQLPNMAWDTIVEQSTANPTLLDSETVKIIANI
IKTNVAVCTSMGADFYPQLGHIYYNMLQLYRAVSSMISAQVAAEGLIATKTPKV
RGLRTIKKEILKLVEYISKARNLDDVVKVLVEPLLNAVLEDYMNNVPDARDAE
VLNCMTTVVEKVGHMIPQGVILILQSVFECTLDMINKDFTEYPEHRVEFYKLLK
VINEKSFAAFLELPAAFKLFVDAICWAFKHNNRDVEVNGLQIALDLVKNIERM
GNVPFANEFHKNYFFIFVSETFFVLTDSDHKSGFSKQALLMKLISLVYDNKIS
VPLYQEAEVPQGTSNQVYLSQYLANMLSNAFPHLTSEQIASFLSALTQYKDLV
VFKGTLRDFLVQIKEVGGDPTDYLFAEDKENALMEQNRLEREKAAKIGGLLKPS
ELDD

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Cytoplasm.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.1

Predicted Molecular Mass: 130.5kDa

Accurate Molecular Mass: 131kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

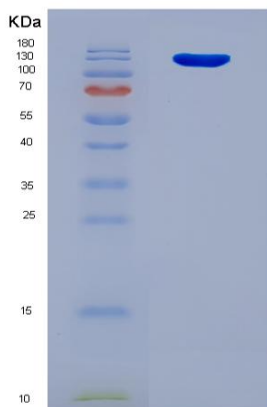


Figure 1. SDS-PAGE

[IMPORTANT NOTE]



The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.