

Recombinant Cytokeratin 18 (CK18)

Catalog No.: **TP10161**

100µg

Sequence Information

Species: Rat

Gene ID: 294853

Swiss Prot: Q5BJY9

Synonyms: KRT18; CYK18; K18; Keratin 18;
Cell proliferation-inducing gene 46
protein; Keratin, type I cytoskeletal 18

Residues: Glu72-Gly380

EKETMQDLNDRLASYLDKVKNLETENRRLESKIREYLEKRGPGQVRDWGHYFKT
IEDLRAQIFANSVDNARIVLQIDNARLAADDFRVKYETELAMRQSVESDIHGLR
KVVDNTNITRLQLETEIEALKEELLFMKKNHHEEVQGLEAQIASSGLTVEVDAP
KSQDLSKIMADIRAQYEQLAQKNREELDKYWSQQIEESTTVVTTKSAEIRDAET
TLLELRRTLQTLTLEIDLDSMKNQINLENLGEVEARYRVQMEQLNGVLLHLESE
LAQTRAEGQRQTQEYEALLNIKVKLEAEIATYRRLLEDG

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Nucleus, Cytoplasm.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.2

Predicted Molecular Mass: 39.6kDa

Accurate Molecular Mass: 40kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

CGATCCGAGAAAGGAGCCATGCAGACCTGAACGACCCGCTGGCCAGCTACCTAGACAAAGGTGAAGAACCCTGGAAACCCGAGAACAGGAGACTAGAGAGCAAAATCCGGGAATATCTCGAGAAAGGGGTCCCCAGGGCGTCAAGAGACTGGGGCCACTACTTCAAGACCATCGAAGATCTC

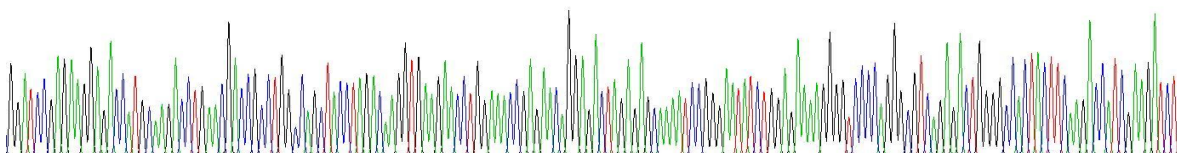


Figure 1. Gene Sequencing (Extract)

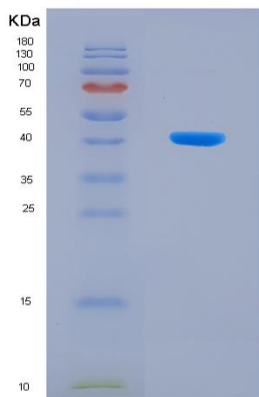


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.