

Recombinant RNA demethylase ALKBH5 (ALKBH5)

Catalog No.: TP10244

100µg

Sequence Information

Species: Human

Gene ID:54890

Swiss Prot:Q6P6C2

Synonyms:ALKB5

Residues:Met1-His394

MAAASGYTDLREKLKSMTSRDNYKAGSREAAAAAAVAAAAAAAAAEPYPVS
GAKRKYQEDSDPERSDYEEQQLQKEEEARKVKSGIRQMRLFSQDECAKIEARID
EVVSRAEKGLYNEHTVDRAPLRNKYFFGEGYTYGAQLQKRGPGQERLYPPGDVD
EIPewHQLVIQKLVEHRVIEGFGVNSAVINDYQPGGCIVSHVDPIHIFERPIV
SVSFFSDSALCFGCKFQFKPIRVSEPVLSLPVRRGSVTVLSGYAADEITHCIRP
QDIKERRAVIILRKTRLDAPRLETksLSSSVLPPSYASDRLSGNNRDPALKPKR
SHRKADPDAAHRPRILEMDKEENRRSVLLPthRRRGsfssenywrksyessedc
SEAAGSPARKVKMRRH

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.4

Predicted Molecular Mass: 45.6kDa

Accurate Molecular Mass: 46kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

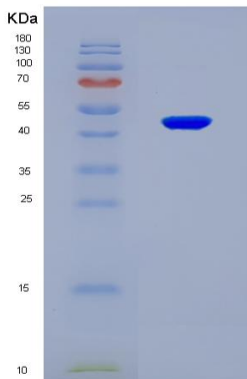


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.