

Recombinant Protein arginine N-methyltransferase 1 (Prmt1)

Catalog No.: TP10275

100µg

Sequence Information

Species: Mouse

Gene ID:15469

Swiss Prot:Q9JIF0-2

Synonyms:ANM1

Residues:Met1-Arg353

MAAAEAANCIMEVSCGQAESSEKPNAEDMTSKDYFDSYAHFGIHEEMLKDEV
TLTYRNSMFHNRHLFKDKVLDVSGTGILCMFAAKAGARKVIGIECSSISDYA
VKIVKANKLDHVVTIIKGKVEEVELPVEKVDIIISEWMGYCLFYESMLNTVLHA
RDKWLAPDGLIFPDRLTYVTAIEDRQYKDYKIHWWENVYGFDMSCIKDVAIKE
PLVDVDPKQLVTNACLIEVDIYTVKVEDLTFTSPFCLQVKRNDYVHALVAYF
NIEFTRCHKRTGFSTSPESPYTHWKQTVFYMEDYLTVKTGEEIFGTIGMRPNAK
NNRDLDFITDLDFKGQLCELSCTDYMR

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Nucleus.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.3

Predicted Molecular Mass: 43.9kDa

Accurate Molecular Mass: 44kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

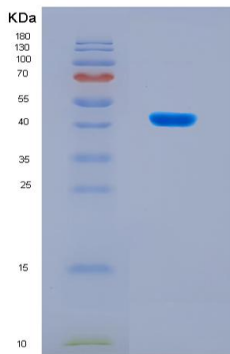


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.