

Recombinant Epoxide hydrolase 1(EPHX1)

Catalog No.: TP11234

100µg

Sequence Information

Species: Human

Gene ID:2052

Swiss Prot:P07099

Synonyms:HYEP

Residues:Met1-Gln455

MWLEILLTSVLGFAIYWFISRDKETLPLEDGGWGPGRSAAREDDSI RPFKVE
TSDEEIHDLHQRIDKFRFTPPLEDS CFHYGFNSNYLKKVISYWRNEFDWKKQVE
ILNRYPHFKTKIEGLDIHFIVHKPPQLPAGHTPKPLLMVHGWPGSFYEFYKIIP
LLTDPKNHGLSDEHVFVVICPSIPGYGFSEASSKKGFNSVATARIFYKLMRLG
FQEFYIQGGDWGSLICTNMAQLVPSHVKG LHLNMAVLVS NFSTLTLLLGQRFRGR
FLGLTERDVELLYPVKEKV FYSLMRESGYMHIQCTKPDTVGSALNDSPVGLAAY
ILEKFSTWTNTEFRYLEDGGLERKFSLDDLLTNVMLYWTTGTI ISSQRFYKENL
GQGWMTQKHERMKVYVPTGFSAFP FELLHTPEKWVRFKYPKLISYSYMVRGGHF
AAFEPELLAQDIRKFLSVLERQ

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag and SUMO

Subcellular Location: Microsome membrane.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.3

Predicted Molecular Mass: 68.9kDa

Accurate Molecular Mass: 69kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

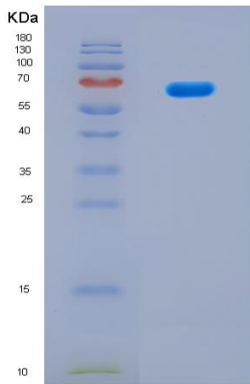


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.