

## Recombinant Serpin A10 (SERPINA10)

Catalog No.: TP10300

100µg

### Sequence Information

**Species:** Mouse

**Gene ID:**217847

**Swiss Prot:**Q8R121

**Synonyms:**PZI; ZPI;

Serine(Cysteine)Proteinase

Inhibitor,Clade A(Alpha-1

Antiproteinase,Antitrypsin)Member 10;

Protein Z-Related Protease Inhibitor;

Protein Z Dependent Protease Inhibitor

**Residues:**Phe22-Leu448

FNLSHSPASVHLESQDYENQTWEYTRTDPREEEEEEEKEEGKDEEYWLRA  
SQQLSNETSSFGFNLLRKISMHDGNVIFSPFGLSVAMVNLMLGKGETKVVQIE  
NGLNLQALSQAGPLILPALFKVKETFSSNRDLGLSQGSFAFIHKDFDIKETYP  
NLSKKYFDIEYVSINFQNSSQARGLINHCIVKETEGKIPKLFDEINPETKLILV  
DYVLFKGGKWLTPFDPSFTEADTFHLDKYRAIKVPMYREGNFTSTFDKKFRCHI  
LKLPHYQGNATMLVVLMEKTGDYLALEDYLTVDLVETWLQNMKTRKMEVFFPKFK  
LNQRYEMHELLKQMGIRRLFSTSADLSELSAMARNLQVSRVLQQSVLEVDERGT  
EAVSGTLSEIIAYSMPPAIKVNRPFHFIIEEMSRMLLFLGRVVNPTVL

### Product Information

**Source:** Recombinant expression.

**Host:** *E.coli*

**Tags:** N-terminal His and GST Tag

**Subcellular Location:** Secreted.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.9

**Predicted Molecular Mass:** 81.8kDa

**Accurate Molecular Mass:** 82kDa as determined by SDS-PAGE reducing conditions.

## **[ USAGE ]**

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

## **[ STORAGE AND STABILITY ]**


**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## **[ IDENTIFICATION ]**

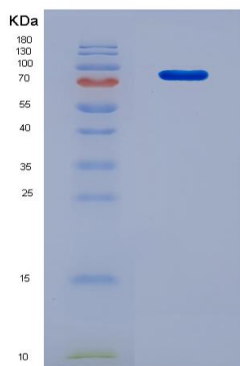

**Figure 1. SDS-PAGE**

## **[ IMPORTANT NOTE ]**

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.