



Recombinant NADH dehydrogenase [ubiquinone] 1 alpha subcomplex subunit 8 (Ndufa8)

Catalog No.: TP11260 100µg

Sequence Information

Species: Mouse

Gene ID:68375

Swiss Prot:Q9DCJ5

Synonyms:NDUA8

Residues:Met1-Val172 W86A

MPGIVELPTLEELKVEEVKVVSAVLKAAAHHYGAQCDKTNKEFMLCRWEKDPK

RCLKEGKLVNGCALNFFRQIKSHCAEPFTEYATCLDYSNMQLFRRHCRQQQAKFD

QCVLDKLGWVRPDLGQLSKVTKVKTDRLPENPYHSRRARPEPNPVIEGDLKPAK

HGTRFFFWTV

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Mitochondrion inner membrane.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.6

Predicted Molecular Mass: 23.2kDa

Accurate Molecular Mass: 23kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

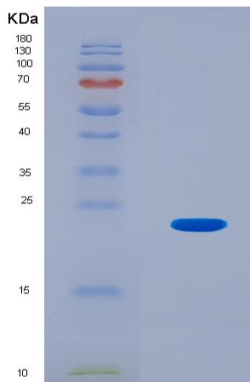


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.