



## Recombinant Human Beta Catenin ( $\beta$ -catenin)

Catalog No.: TP10320

50 $\mu$ g

### Sequence Information

**Species:** Human

**Gene ID:**1499

**Swiss Prot:**P35222

**Synonyms:**CTNNB1

**Residues:**Ala2-Val451

ATQADLMELDMAMEPDRKAAVSHWQQSYLDSGIHSGATTTAPSLSGKGN  
PEEEDVDTSQVLVEWEQGFSSQFTQEQVADIDGQYAMTRAQRVRAAMFPE  
TLDEGMQIPSTQFDAAHPTNVQRLAEPQMLKHAVVNLINYQDDAELATR  
AIPELTKLLNDEDQVVVNKAAMVHQLSKKEASRHAIMRSPQMVSIVRT  
MQNTNDVETARCTAGTLHNLSSHREGLLAIFKSGGIPALVKMLGSPVDSV  
LFYAITTLHNLHLEHQAQMAVRLAGGLQKMAVLLNKTNVKFLAITTDCL  
QILAYGNQESKLIIASGGPQALVNIMRTYTYEKLWTTSRVLKVLVSVCS  
SNKPAIVEAGGMQALGLHLDPSQRLVQNCLWTLRNLSDAATKQEGMEGL  
LGTLVQLLGSDDINVVTCAAGILSNLTCNNYKNKMMVCQVGGIEALVRTV

### Product Information

**Source:** Prokaryotic expression.

**Host:** *E. coli*

**Tags:** N-terminal His and GST Tag

**Subcellular Location:** Membrane, Nucleus, Cytoplasm.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS (PH7.4) , containing 5% Trehalose.

**Original Concentration:** 200 $\mu$ g/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.9

**Predicted Molecular Mass:** 79.1kDa

**Accurate Molecular Mass:** 79kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in PBS (pH7.4) to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.


Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The



loss rate is less than 5% within the expiration date under appropriate storage condition.

## **[ IDENTIFICATION ]**

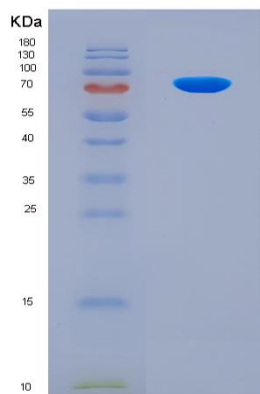

**Figure . SDS-PAGE**

## **[ IMPORTANT NOTE ]**

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.