

Recombinant Lymphotoxin Beta (LTb)

Catalog No.: TP11268

100µg

Sequence Information

Species: Mouse

Gene ID:16994

Swiss Prot:P41155

Synonyms:LT-B; p33; TNFC; TNF-C; TNFSF3

Residues:Ile58-Met304

IIIGSGAQAQKRLDDSKPSCILPSPSSLSETPDRLHPQRSNASRNLASTSQGPVAQSSREASAWMTILSPAADSTPDGPVQQLP
KGEPETDLNPELPAAHLIGAWMSGQGLSWEASQEEAFLRSGAQFSPTHGLALPQDGVIYLYCHVGYRGRTPPAGRSRARSLLTR
SALYRAGGAYGRGSPELLLEGAETVTPVVDPIGYGSLWYTSVGFGLAQLRSGERVVNIHPDMVDYRRGKTFFGAVM

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: Two N-terminal Tags

Subcellular Location: Membrane, Chromosome.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.0

Predicted Molecular Mass: 28.4kDa

Accurate Molecular Mass: 32kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is

less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

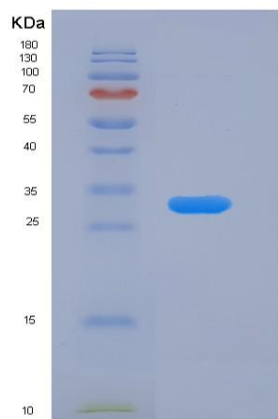


Figure . SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.