

Recombinant Cluster Of Differentiation 99 (CD99)

Catalog No.: **TP11273**

50µg

Sequence Information

Species: Human

Gene ID:4267

Swiss Prot:P14209

Synonyms:MIC2X; MIC2Y; MIC2; Antigen

Identified By Monoclonal Antibodies

12E7,F21 And O13; E2 antigen; T-cell
surface glycoprotein E2

Residues:Asp23~Asp122

DGGFDLSDALPDNENKKPTAIPKKPSAGDDFDLGDAVVDGENDDPRPPNPPKPM

PNPNNHPSSSGSFSDADLADGVSGGEGKGGSDGGGSHRKEGEEAD

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: C-terminal His and Fc Tag

Subcellular Location: Membrane.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 60µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.0

Predicted Molecular Mass: 26.6kDa

Accurate Molecular Mass: 50kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

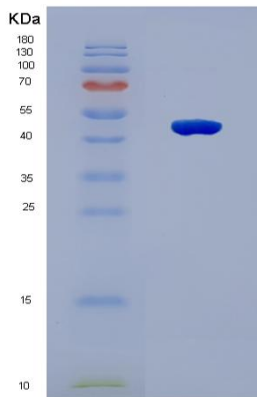


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.