



Recombinant Human Tissue-Type Plasminogen Activator/PLAT Protein

Catalog No.: TP10329 100µg

Sequence Information

Species: Human

Gene ID:5327

Swiss Prot:P00750

Synonyms: T-PA; TPA; t-plasminogen activator; Tissue plasminogen activator

Residues:Ser36-Pro562

SYQVICRDEKTQMIYQQHQSWLRPVLRSNRVEYCWCNSGRAQCHSVPVKSCSEPRCFNGGTCQQALYFSDFVCQCPEGFAGKC
CEIDTRATCYEDQGISYRGTWSTAESGAECTNWNSSALAQKPYSGRRPDAILRLGLGNHNYCRNPDRDSKPWCYVFKAGKYSSE
FCSTPACSEGNSDCYFGNGSAYRGTHSLTESGASCLPWNSMILIGKVYTAQNPSAQALGLGKHNYCRNPDGDAKPWCHVLKNR
RLTWEYCDVPSCSTCGLRQYSQPQFRIKGGFLFADIASHPWQAAIFAKHRRSPGERFLCGGILISSCWILSAAHCFQERFPPHH
LTVILGRITYRVVPGEEEEQKFEVEKYIVHKEFDDDTYDNDIALQLKSDSSRCAQESSVVRTVCLPPADLQLPDWTECELSGYG
KHEALSPFYSERLKEAHVRLYPSSRCTSQHLLNRTVTDNMLCAGDTRSGGPQANLHDACQGDSSGGPLVCLNDGRMTLVGIISW
GLGCGQKDVPGVYTKVTNYLDWIRDNMRP

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: C-terminal His-Tag

Subcellular Location: Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS (PH7.4) , containing 5% Trehalose.

Original Concentration: 400µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.1

Predicted Molecular Mass: 60.1kDa

Accurate Molecular Mass: 62kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in PBS (pH7.4) to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

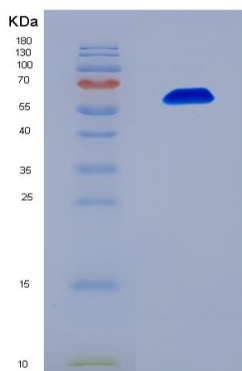


Figure . SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.