



Recombinant LINE-1 retrotransposable element ORF2 protein (LORF2)

Catalog No.: TP11235 100µg

Sequence Information

Species: Human

Gene ID: O00370

Swiss Prot: O00370

Synonyms: ORF2p

Residues: Met1-Ile239

MTGSNSHITILT LNVNGLNSPIKRHRLASWIKSQDPSVCCIQTHTLCRDTHRLKIKGWRKIYQANGKQKKAGVAILVSDKTDFKPTKIKRDKEGHYI
MVKGSIQQEELTILNIYAPNTGAPRFIKQVLSDLQRDLDSHTLIMGDFNTPLSILDRSTRQKVNKDTQELNSALHQTDLIDIYRTLHPKSTEYTFSSA
PHHTYSKIDHIVGSKALLSKCKRTEIITNYLSDHSAIKLELRI

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag+GST

Subcellular Location: /.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose .

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 10.2

Predicted Molecular Mass: 55.6kDa

Accurate Molecular Mass: 58kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

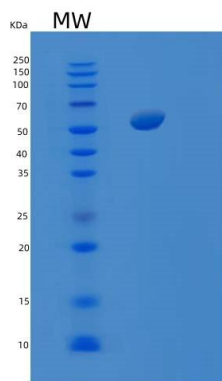


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.