

Recombinant Human Reversion-inducing cysteine-rich protein with Kazal motifs(RECK)

Catalog No.: TP11327

50µg

Sequence Information

Species: Human

Gene ID:8434

Swiss Prot:O95980

Synonyms:

Residues:Gly27-Pro941

GLAPGSAGALCCNHSKDNQMCRDVCEQIFSSKSESRLKHL LQRAPDYCPETMVE
IWNCMNSSLPGVFKKSDGWGLGCC ELAIALECRQACKQASSKNDISKVCRKEY
ENALFSCISRNEGMSVCCSYAGHHTNCREYCQAI FRTDSSPGPSQIKAVENYCA
SISPLIHCVNNTQSYPMRNP TDSLYCCDRAEDHACQNACKRI LMSKKT EMEI
VDGLIEGCKTQPLPQDPLWQC FLESSQSVHPGVTVHPP PSTGLDGAKLHCCSKA
NTSTCRELCTKLYSMWGNTQSWQEFDRFCEYNPVEVSMLTCLADVREPCQLGC
RNLTYCTNFNNRPT ELFRSCNAQSDQGAMNDMKLWEKGS IKMPF INIPVLDIKK
CQPEMWKAIACSLQIKPCHSKSRGSIICKSDCVEILKKCGDQNKFPEDHTAESI
CELLSPTDDLKNCIPLD TYLRPSTLGNIVEEVTHPCNPNPCAN ECVNRKGC
PSGDPLPYFCVQGCKLGEASDFIVRQGT LIQVPSSAGEVG CYKICSCGQSGLL
ENCMEMHCIDLQKSCIVGGKRKSHGTSFSIDCNVCSCFAGNLVCSTR LCLSEHS
SEDDRRTFTGLPCNCADQFVPVCGQNGRTYPSACIARCVGLQDHQFEFGSCMSK
DPCNPNPCQKNQRCIPKPQVCLTTFDKFGCSQYECVPRQLACDQVQDPVCDTDH
MEHNNLCTLYQRGKSLSYKGPCQPF CRATEPVCGHNGETYSSVCAAYS DRVAVD
YYGDCQAVGVLSEHSSVAECASVKCPSLLAAGCKPIIPP GACCPLCAGMLRVLF
DKEKLDTIKVTNKKPITVLEILQKIRMHVSV PQCDVFGYFSIESEIVILIIPV
DHYPKALQIEACNKEAEKIESLINS DSPTLASHVPLSALIISQVQVSSSV

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: C-terminal His Tag

Subcellular Location: Membrane.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000 μ g/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.5

Predicted Molecular Mass: 102kDa

Accurate Molecular Mass: 103kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

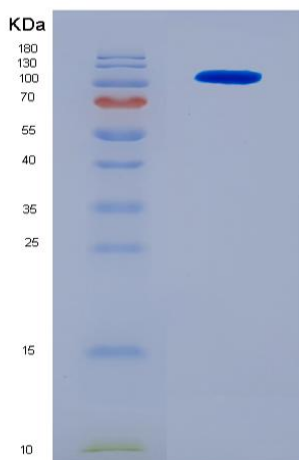


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.

