



Histone H2A.X Recombinant Rabbit mAb

RMA02701

100ul

[PROPERTIES]

Source: Monoclonal antibody preparation

Host: Rabbit

Purification: Protein A affinity chromatography.

Traits: Liquid

Concentration: 0.8mg/ml

Clone Number:C1

Ig Isotype:IgG

UOM: 100ul

Cross Reactivity: Human, Mouse, Rat

Applications: WB, IHC, ICC/IF, FC, IP

[IMMUNOGEN]

Immunogen: A synthetic peptide of human Histone H2A.X

[APPLICATIONS]

WB: 1:1000-1:5000

IHC: 1:200-1:1000

ICC/IF: 1:200-1:1000

FC: 1:20

IP: 1:20-1:100

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.




Store at 4°C for frequent use.

Aliquot and store at -20°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

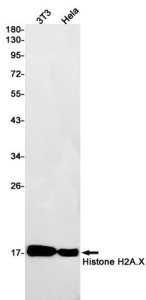

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.