

Recombinant FK506 Binding Protein 10 (FKBP10)

Catalog No.: TP11507

100µg

Sequence Information

Species: Human

Gene ID:60681

Swiss Prot:Q96AY3

Synonyms:FKBP65; hFKBP65; Rotamase; 65

kDa FK506-binding protein;

Immunophilin FKBP65; Peptidyl-prolyl

cis-trans isomerase FKBP10

Residues:Val158-Asp440

DTVQVSTLLRPPHCPRMVQDGFVRYHYNGTLLDGTSEFDTSSKGGTYDTYVGS

GWLIKGMDDQGLLGMPGERRKIIIPPFLAYGEKGYGTVIPQASLVFHVLLIDV

HNPKDAVQLETLLEPPGCVRRAGAGDFMRYHYNGSLMDGTLFDSSYSRNHTYNT

YIGQGYIIPGMDQGLQGACMGERRRITIPPHLAYGENGTGDKIPGSAVLIFNVH

VIDFHNPAADVVEIRTLRPSSETCNETTKLGDFVRYHYNCSSLDDGTQLFTSHDYG

APQEATLGANKVIEGLD

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His Tag

Subcellular Location: Endoplasmic reticulum lumen .

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.1

Predicted Molecular Mass: 32.5kDa

Accurate Molecular Mass: 33kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

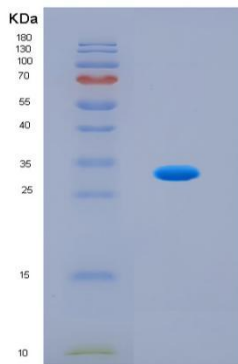


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.