

## Recombinant Human HSP70 Protein

Catalog No.: TP02525

100µg

### Sequence Information

**Species:** Human

**Gene ID:**3303/3304

**Swiss Prot:**P0DMV8

**Synonyms:**Heat shock 70 kDa protein 1,  
HSP70.1, HSP70-1/HSP70-2, HSPA1A,  
HSPA1, HSPA1B

**Residues:**Met1-Asp641

MAKAAAIGIDLGTTSYSCVGVFQHGKVEIANDQGNRTTPSYVAFTDTERLIGDA  
AKNQVALNPQNTVFDARLIGRKFQDPVQSDMKHWPQVINDGDKPKVQVSYK  
GETKAFYPEEISSMVLTKMKEIAEAYLGYPVTNAVITVPAYFNDSQRQATKDAG  
VIAGLNVLRIINEPTAAAIAYGLDRTGKGERNVLIFDLGGGTFDVSILTIDGDI  
FEVKATAGDTHLGGEDFDNRLVNHVVEEFKRKHKKDISQNKRAVRRRLRTACERA  
KRTLSSTQASLEIDSLFEGIDFYTSITRARFEELCSDLFRSTLEPVEKALRDA  
KLDKAQIHDLVLVGGSTRIPKQKLLQDFNGRDLNKSINPDEAVAYGAAVQAA  
ILMGDKSENVQDLLLLDVAPLSLGLLETAGGVM TALIKRNSTIPTKQTQIFTTYS  
DNQPGVLIQVYEGERAMTKDNNLLGRFELSGIPPAPRGVQIEVTFDIDANGIL  
NVTATDKSTGKANKITITNDKGRLSKEEIERMVQEAKEYKAEDVQRERVSANK  
ALESYAFNMKSAVEDEGLGKGISEADKKKVLDKCQEVISWLDANTLAEKDEFEH  
KRKELEQVCNPIISGLYQGAGGPGGGFGAQGPKGGSGSGPTIEEVD

### Product Information

**Source:** Recombinant expression.

**Host:** *E.coli*

**Tags:** N-terminal His-Tag

**Subcellular Location:** Cytoplasm.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.3

**Predicted Molecular Mass:** 72.2kDa

**Accurate Molecular Mass:** 72kDa as determined by SDS-PAGE reducing conditions.

## [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

## [ STORAGE AND STABILITY ]


**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

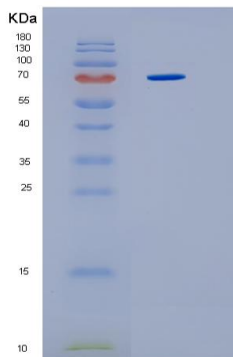

**Figure 1. SDS-PAGE**

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.