

Recombinant Human JNK1 / MAPK8 Protein

Catalog No.: TP06466

50µg

Sequence Information

Species: Human

Gene ID:5599

Swiss Prot:P45983

Synonyms:JNK;JNK-46;JNK1;JNK1A2;
JNK21B1/2;PRKM8;SAPK1;SAPK1c

Residues:Met1-Arg427

MSRSKRDNNFYSVEIGDSTFTVLKRYQNLKPIGSGAQGIVCAAYDAILERNVAI
KKLSRPFQNQTHAKRAYRELVLMKCVNHKNIIGLLNVFTPQKSLEEFQDVYIVM
ELMDANLCQVIQMELDHERMSYLLYQMLCGIKHLHSAGIIHRDLKPSNIVVKS
CTLKILDFGLARTAGTSFMMPYVWTRYRAPEVILGMGYKENVDLWSVGCIMG
EMVCHKILFPGRDYIDQWNKVIEQLGTPCEFMKKLQPTVRTYVENRPKYAGYS
FEKLFDPVLFPADEHNKLLKASQARDLLSKMLVIDASKRISVDEALQHPYINW
YDPSEAEAPPPKIPDKQLDEREHTIEEWKELIYKEVMDLEERTKNGVIRGQSP
LGAAVINGSQHPSSSSSVNDVSSMSTDPTLASDTSLEAAAGPLGCCR

Product Information

Source: Eukaryotic expression.

Host: *Baculovirus-Insect Cells*

Tags: N-terminal GST-Tag

Subcellular Location: Cytoplasm.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.9

Predicted Molecular Mass: 75.0kDa

Accurate Molecular Mass: 75kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

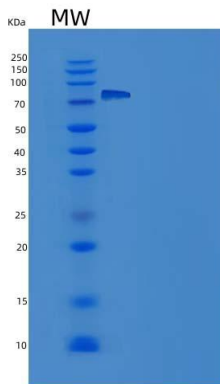


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.