



Notch1 Recombinant Rabbit mAb

RMA01263 100ul

[PROPERTIES]

Source: Monoclonal antibody preparation

Host: Rabbit

Purification: Protein A and Protein G affinity chromatography.

Traits: Liquid

Concentration: 0.5mg/ml

Ig Isotype:IgG

UOM: 100µl

Cross Reactivity: Human, Mouse

Applications: WB,IHC,ICC/IF,FC.

[IMMUNOGEN]

Immunogen: Recombinant protein of human Notch1.

[APPLICATIONS]

WB: 1:1000

IHC-P: 1:100-1:200

ICC/IF: 1:100-1:200

FC: 1:100-1:200

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.



Store at 4°C for frequent use.

Aliquot and store at -20°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

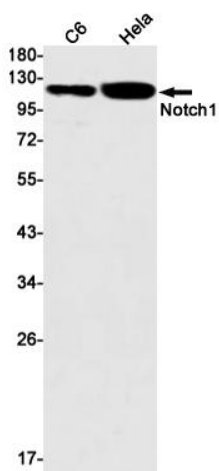


Figure. Western Blot; Sample:
Recombinant Notch1, Human.

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.