

Recombinant Caspase 3 (CASP3)

Catalog No.: TP11743

50µg

Sequence Information

Species: Human

Gene ID:836

Swiss Prot:P42574

Synonyms:CPP32; CPP32B; SCA1; Apoptain;

Yama; Apoptosis-Related Cysteine

Peptidase; CysteinyI Aspartate Specific

Proteinases 3; SREBP cleavage activity 1

Residues:Ser176-His277

SGVDDDMACHKIPVEADFLYAYSTAPGYYSWRNSKDGSWFIQSLCAMLKQYADK

LEFMHILTRVNRKVATEFESFSFDFATFHAKKQIPCIIVSMLTKELYFYH

Product Information

Source: Prokaryotic expression.

Host: *E. coli*

Tags: No Tag

Subcellular Location: Cytoplasm.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: 20mM Tris, 150mM NaCl, pH8.0, containing 0.01% skl,5%Trehalose.

Original Concentration: 50µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.9

Predicted Molecular Mass: 11.9kDa

Accurate Molecular Mass: 12 kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was

determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

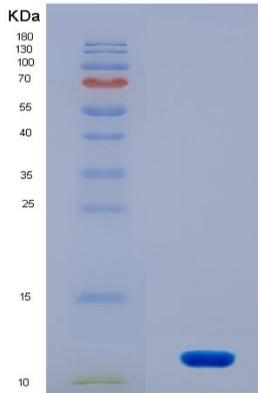


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.