

## Recombinant Angiotensin Like Protein 7 (ANGPTL7)

Catalog No.: TP10842

100µg

### Sequence Information

**Species:** Mouse

**Gene ID:**654812

**Swiss Prot:**Q8R1Q3

**Synonyms:**CDT6; AngX; Cornea-Derived  
Transcript 6; Angiotensin-related  
protein 7; Angiotensin-like factor

**Residues:**Gln22-Pro337

QKPHKRKTQLKAAGCCEEMRELKAQVANLSSLLGELSRKQESDWVSVMQVMEL  
ESSSKHMESRLSTAESKYSEMNNQIDIMQLQAAQTVTQTSADAIYDCSSLYQKN  
YRISGVYKLPDEF LGSPELEVFCDMETSGGGWTIIQRRKSGLVSFYQDWRQYK  
QGFGSIRGDFWLNEHHRHLTRQPSRLRVELEDWEGNARYAEYSYFALGNELNS  
YRLF LGNYSGNVGDALLYHNNTVFSTKDKDNDNCLDKCAQLRKGGYWYNCCTD  
SNLNGVYYRLGEHRKHMDGISWYGWHGANYS LKR VEMKIRPEAFKP

### Product Information

**Source:** Recombinant expression.

**Host:** *E.coli*

**Tags:** N-terminal His Tag

**Subcellular Location:** Secreted.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 7.9

**Predicted Molecular Mass:** 39.9kDa

**Accurate Molecular Mass:** 37kDa as determined by SDS-PAGE reducing conditions.

### **[ USAGE ]**

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### **[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

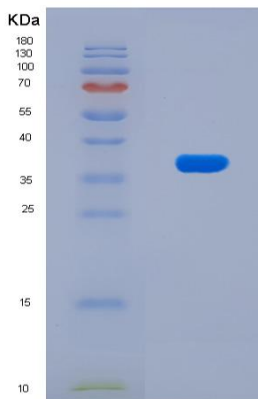


Figure 1. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.